

CHEMTRIX
Scalable Flow Chemistry

**Flow Chemistry
Method Development**



GramFlow[®]

Flow Chemistry Method Development

GramFlow[®] is a glass continuous flow reactor designed for laboratory based screening & reaction optimization. With its integrated pre-heat & heat exchange channels, GramFlow[®] has optimal thermal control for challenging flow chemistry applications of the A+B type. The reactor is suited to a wide range of chemical applications;

- Performance of reaction screening & optimization
- Assessment of process feasibility
- Perfect entry-level research reactor
- Material production at the g-scale

GLASS REACTORS

- Excellent heat & mass transfer
- Excellent mixing using zig-zag structure
- German quality

SPECIFICATIONS


- Reaction types: A + B → P
- Throughput: 0.2 to 10 ml/min (up to 600 g/h)
- Reaction volume: 1 ml
- High operating pressure: 20 bar
- Wide temperature range: -20 to 150°C
- Inert materials: PTFE, FFKM, Glass

DIMENSIONS

- 126 (W) x 61 (D) x 46 mm (H)



CHEMTRIX BV DEVELOPS & PRODUCES
CONTINUOUS FLOW REACTORS & SYSTEMS
FROM LAB TO PRODUCTION

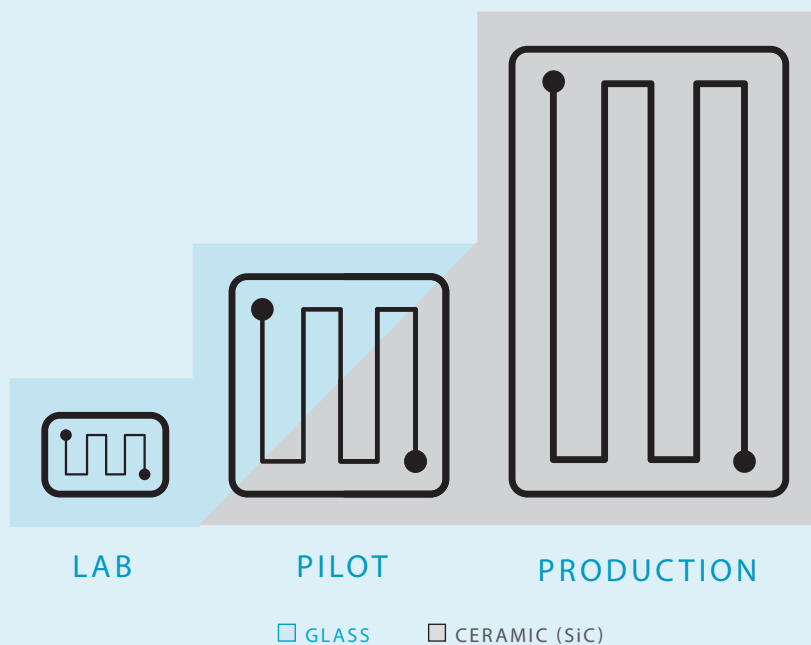



CHEMTRIX BV IS HEADQUARTERED IN THE NETHERLANDS

WITH OUR GLOBAL OFFICES & LABORATORIES
WE ASSIST OUR CUSTOMERS WITH LOCAL CHEMICAL
& TECHNICAL SUPPORT

CHEMTRIX
Scalable Flow Chemistry

Chemtrix BV
Galvaniweg 8a
6101 XH Echt
The Netherlands

T +31 46 70 22 600
info@chemtrix.com
www.chemtrix.com

strategic partner of
3M Science.
Applied to Life.™

Details of our global offices and labs can be found at www.chemtrix.com/contact

The data in this brochure was valid at the time of publishing, and is subject to change without notice. Chemtrix BV reserves the right to modify its product at any time without further obligation.